



„Buttockmate“ type 2 - Intramuscular Injection Simulator

Order code: **4109.M153-1**

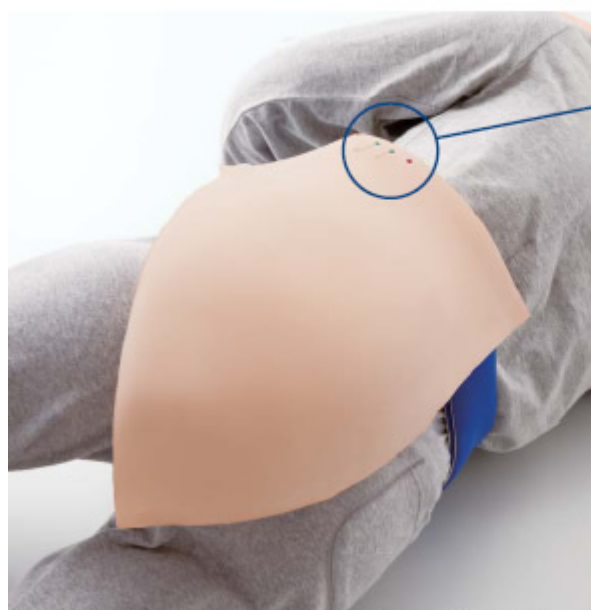

The epidermis is made of a special material that leaves no needle mark.

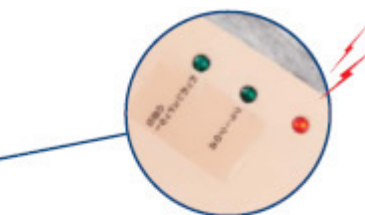

Selection of an incorrect site is indicated by blinking of the red lamp and buzzing. Successful puncturing at the correct site is indicated by the green lamp.
 (Only at Clarke's point and Hochstetter's point)



Being compact, the device needs a small space for storage.

Cena bez DPH

2.435,00 Eur

Price with VAT

2.946,35 Eur

Parameters

Cannula, injection, puncturing - filter

Buttock

Quantitative unit

ks

Reproducing Clarke's point and Hochstetter's point. Allowing simulation in the lateral recumbent position.

Skeletal reference points such as the anterior superior and anterior posterior iliac spines and greater trochanter can be felt for efficient training in finding the injection site. Fluid can be injected. Do not use a drug solution. Use water (tap water only) for practice.

- **Weight:** about 1,2 kg
- **Accessories:** Belt, Drain bag, Storage case